

ULTRASONIC CLEANERS WITH SWEEP TECHNOLOGY

Quala Ultrasonic Cleaners feature self-contained generators in a variety of sizes and capacities to provide efficiency and dependability for all your cleaning applications. Constructed with a durable stainless steel tank, with vinyl-clad steel and aluminum housing. Equipped with a 60-minute timer. 2½ year warranty on all cleaners. Quala Sweep Machines utilize sweep technology, a cleaning wave that sweeps back and forth throughout the ultrasonic tank, reaching all corners of the machine - no reduced cleaning areas.

FEATURES

Swivel drain

Cover included

BENEFITS

Convenient system for draining the machine

Reduces the possibility of water splashing to work area

CODE	DESCRIPTION
------	-------------

5200 Sweep with Timer

1209	Internal tank dimensions: 11¾"x 6"x 4" (29.8 x 15.2 x 10.2 cm) Capacity: 4¼ quarts (4 liters)
-------------	--

5200 Accessory Kit

1212	Includes one full size stainless steel basket, positioning cover, two 600 mL beakers complete, one bur tray and one bur box
-------------	---

5300 Sweep with Timer

1210	Internal tank dimensions: 15½"x 9"x 6" (39.4 x 22.9 x 15.2 cm) Capacity: 13 quarts (12.3 liters)
-------------	---

5300R Sweep (Recessed) with Timer

1211	Internal tank dimensions: 15½"x 9"x 6" (39.4 x 22.9 x 15.2 cm) Capacity: 13 quarts (12.3 liters)
-------------	---

5300 & 5300R Accessory Kit

1213	Includes one solid side basket, positioning cover, two 600 mL beakers complete, one bur tray and one bur box
-------------	--



5200




5300




5300R

SWEEPZONE Performance Power

Automatically controlled frequency changes create an ultrasonic cleaning "grid" which literally "sweeps" through the entire tank, unlike that pattern created by traditional ultrasonic systems. The result is total cleaning, which means constant power and consistent output.

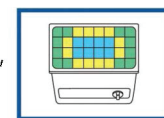

Traditional Ultrasonic System

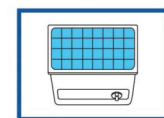

SweepZone Ultrasonic System

- Diminished cleaning power
- Lowest cleaning power
- Indicates areas of constant, intense cleaning power

Distributed by:



Quala
DENTAL PRODUCTS
www.quala.com

#QSS074R2